


Never Tear Another Sinus Membrane

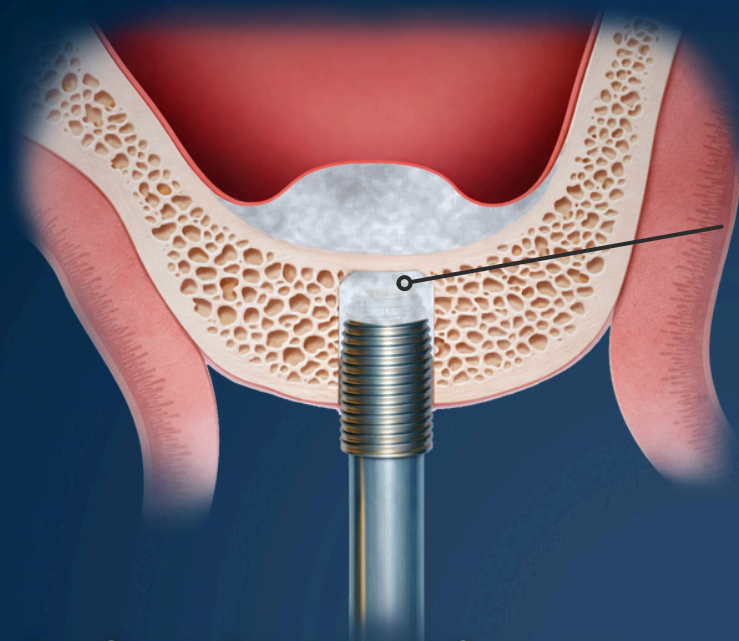

BioDensification is designed to flow through porous cancellous bone channels and lift the sinus membrane

For Sinus Bumps and Complete Sinus Augmentation

TOSA™

Trans-Osseous Sinus Augmentation™

- Stop your pilot drill short of penetrating into the sinus
- Prepare the osteotomy with the **TOSA™ Drill** short of the sinus
- Fill osteotomy with **BioDensification™** and use the driver
- The TOSA™ Driver acts as a piston, generating controlled pressure that transmits through trabecular spaces to gently elevate the sinus membrane
- After lifting membrane, penetrate the sinus to add additional bone graft if needed
- Place the implant at the time of grafting



TOSA™ Kit:

Drill and Driver (patent pending)
BioDensification sold separately

STEINERBIO®

FDA-cleared bone grafts.
Over a decade of clinical use.

Learn more at steinerbio.com/TOSA

